# *aquapick*<sup>®</sup>

CORDLESS ORAL IRRIGATOR

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Advanced Oral Irrigation System  
Sistem avansat de irigare orală

**AQ-208**

# IMPORTANT SAFE GUARDS

## 1. DANGERS

- ① If you do not know the operating process and safety instructions, please do not rush to use the product to avoid hurts caused by mis-operation.
- ② Do not disassemble the product or change the parts of the product to avoid damage.
- ③ Do not use water above the temperature of 50°C.
- ④ If the USB cable is damaged, please stop using and contact service center.
- ⑤ Do not put the product in a fire or near a heat.
- ⑥ Also do not charge, use, or leave it at high temperature.

## 2. WARNINGS

- ① Do not use this product for any other purpose other than oral cleaning.
- ② Do not push the jet tip hard on the gums or teeth to avoid hurting on the gums.
- ③ Do not immerse the product in water for a long time; do not cause fire, electric shock to avoid operation failure.
- ④ Do not damage, destroy or modify the USB cable, and do not stretch, bend or twist the USB cable. Do not press heavy objects on the USB cable or get it caught between objects, and keep the product and cable away from hot objects.
- ⑤ Do not store the product into a hot or humid places, and avoid direct sunlight.
- ⑥ Do not charge the product in hot or humid place, such as the bathroom or toilet.
- ⑦ If not in use of the product for a long time, please unplug the USB cable to avoid fire or danger.

## 3. NOTICE

- ① Do not use without a jet tip.
- ② This product is not suitable for children and infants unless they are supervised or guided by their safety responsible supervisor on how to use this product. Keep out of reach for children or infants. Do not use it as a toy.
- ③ If you have periodontal disease or have had oral surgery in the past 2 months, please consult your dentist before using.
- ④ Do not add mouthwash into the water tank. This will cause corrosion and the water tank will be adhered to the body resulting in the water tank being stuck.

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

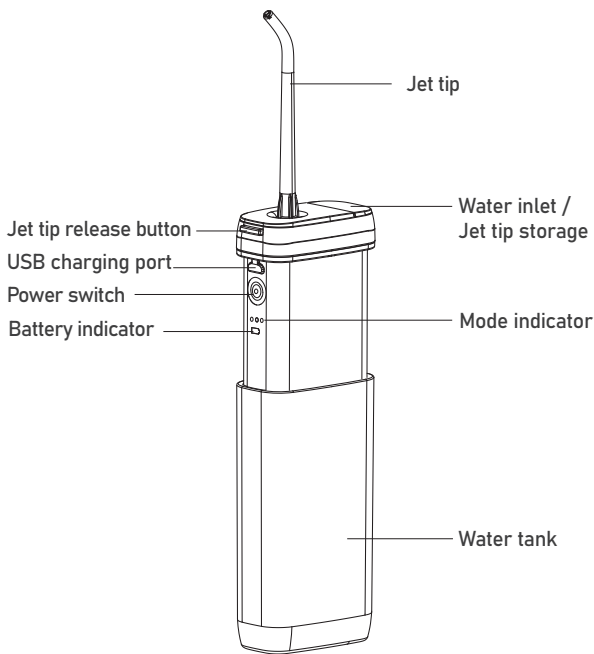
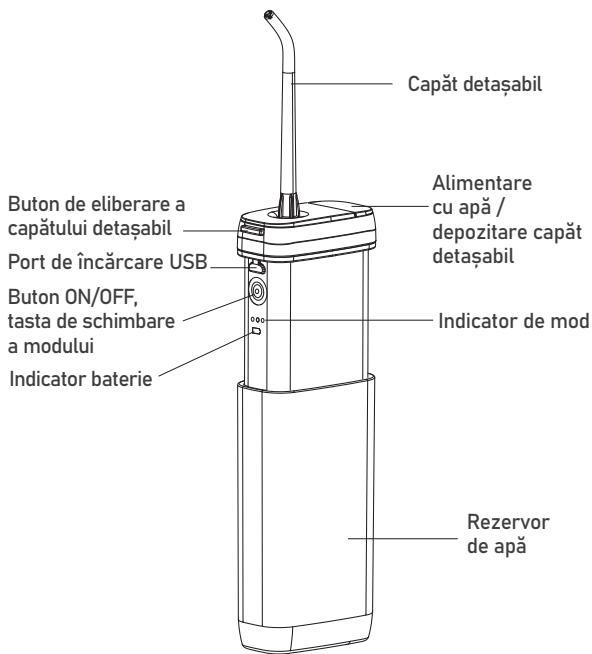
Install the jet tip.

2

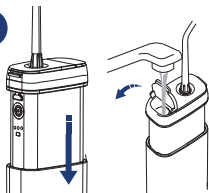Pull out the water tank from the body, open the water inlet, fill it with water, and then close it.

3

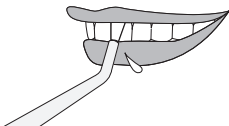Put the jet tip in your mouth and press the power button to start irrigating.

4


Long press the power switch button to switch the mode.

5

For storage after use, remove the jet tip, open the water inlet, and put the jet tip inside.

6

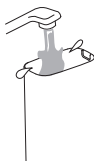Regular cleaning : detach the water tank and clean the inner part.

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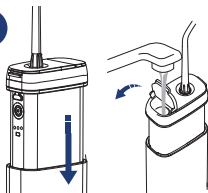
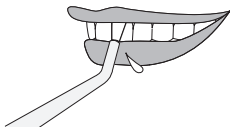
## HOW TO CHANGE MODE

Mode	Mode Switch	Mode indicator color	Target user	Water pressure
Soft Mode	Press the switch to start flossing	Green	Gentle water pressure, gentle cleaning, suitable for sensitive teeth and first time users	40~60 PSI
Massage Mode	Press the switch and hold for 0.5' to switch to pulse mode	Blue	Variable frequency pulse pressure alternate the strong and the weak mode can deeply clean and massage the gums	60~100 PSI
Normal Mode	Press the switch and hold for 0.5' to switch to pulse mode	White	Strong water pressure, efficient cleaning, meeting the daily cleaning for the most of the people	120~140 PSI
Battery indicator	<ol style="list-style-type: none"> <li>1. The indicator turns red when charging is insufficient.</li> <li>2. The indicator flashes red when charging.</li> <li>3. When fully charged, the indicator turns green.</li> </ol>			

※ For more information, contact the Aquapick Customer Service Center.

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## PRODUCT MAINTENANCE

### 1. Jet tip

Press the jet tip button and remove the jet tip, wash it with water and then wipe dry with a cloth.

Regular jet tips should be replaced every 3 months. If the top of the jet tip is deformed, replace it with a new one to avoid malfunction.

### 2. Main body

Remove the dirt with a cloth dampened with liquid soap, then wipe clean with a dry cloth. If the product body is dirty, wash it with water.

Do not wipe the product with thinner, gasoline, alcohol etc. Do not wash with hot water above 50°C.


Please wipe dry after cleaning.

### 3. Water inlet

Wipe off the dirt, otherwise it will cause a leakage.

### 3. Water tank

Open the tank and clean it regularly  
Wash the internal with a brush or other object and wipe with a dry cloth  
When not in use of the product for a long time, please keep the water tank dry. Thank you for purchasing our product.

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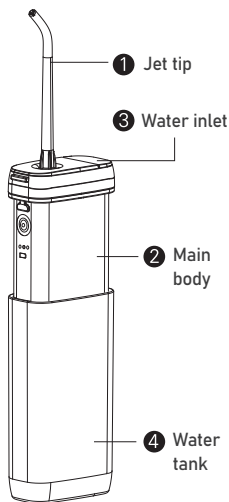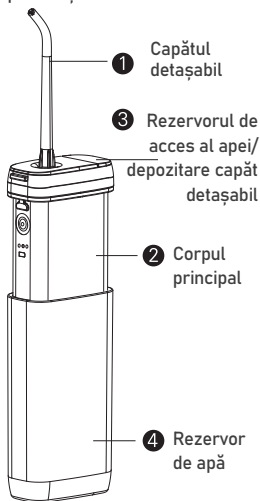
## TROUBLESHOOTING

Problem	Possible cause	Solution
Unit does not work	Newly purchased or not used for more than three months	Please charge it for more than 4 hours
Low water pressure when working	Something is stuck inside the jet tip	Replace the jet tip
Battery Only can be used for a short time when fully charged	The battery is at end of its life cycle	Contact your seller
	Insufficient charging time	Please charge it for more than 4 hours before use
The battery indicator does not light up	The power plug is not properly inserted.	Please connect the USB cable correctly.
	Product USB cable circuit is faulty	Please contact your seller

### Note

Because the gum is not used to the product's stimulation, the gum may bleed at the initial stage. With continuous use with low water pressure, the bleeding usually stops within 2 weeks. If the bleeding lasts for more than two weeks, please stop using the product and consult your dentist.

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