# *aquapick*<sup>®</sup>

## IRIGATOR BUCAL


### User Guide

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
AQ-320

## PRODUCT STRUCTURE

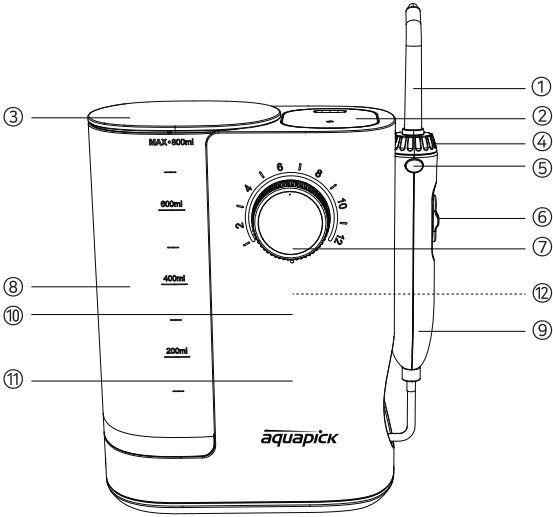- |                          |                            |
|--------------------------|----------------------------|
| ① Jet Tip                | ⑦ Water Pressure Switch    |
| ② Jet Tip Holder Lid     | ⑧ Water Tank               |
| ③ Water Tank Lid         | ⑨ Handle                   |
| ④ Tip Rotating Knob      | ⑩ Jet Tip Holder           |
| ⑤ Jet Tip Release Button | ⑪ Main Body                |
| ⑥ Handle Switch          | ⑫ (Built in) UV sterilizer |

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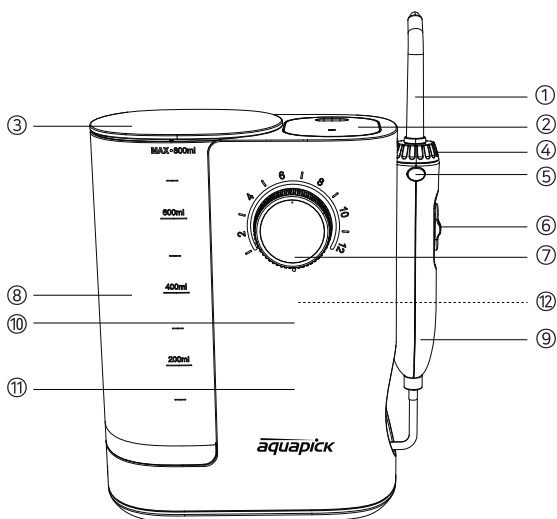
## HOW TO USE

Plug the adaptor into an electrical outlet.

### 1. Storing jet tips and UV sterilization

You can use the jet tip holder in two ways. The first is to use it wherever you want without putting it in the UV sterilizer inside the body. The second is to put it in the UV sterilizer inside the body and sterilize the jet tips. When you close the lid of the jet tip holder, the UV sterilizer turns on, sterilizes jet tips, and turns off automatically after 15 minutes. If you open the jet tip holder lid during sterilization, UV sterilization is automatically stopped to protect users from UV rays, and the UV sterilizer indicator light is also turned off. (Caution: keep the jet tip holder dry)

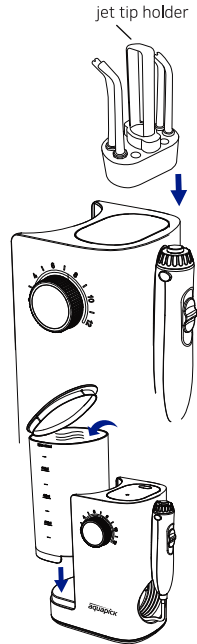
### 2. Filling Water

Remove the water tank and water tank lid from the body and open the water tank lid. Fill water into the water tank and put back the water tank lid to its original place. Then place the water tank back onto the base and press it down gently. Ensure that the rubber seal of the water tank is aligned with the rubber seal of the device base.

### 3. Inserting And Removing The Tip

Insert the jet tip at the hole on the top part of the handle. Press down until there is a clicking sound. To remove the jet tip from the handle, press the jet tip release button on the top part of the handle and gently pull out the jet tip.

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
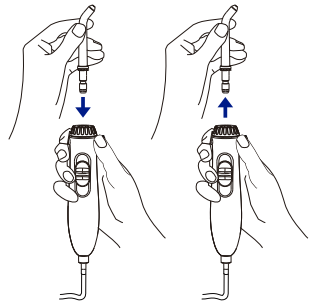
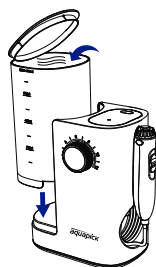
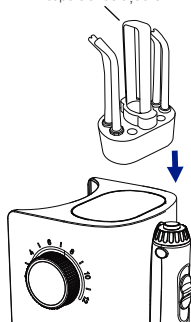
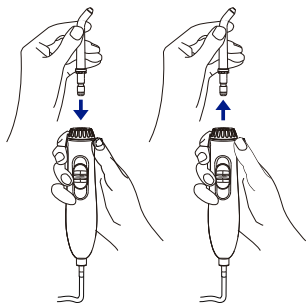
#### 4. Adjusting Water Pressure

The water pressure can be adjusted in 12 levels. Turning the water pressure switch clockwise will increase the water pressure and turning the water pressure switch counter clockwise will decrease the water pressure.
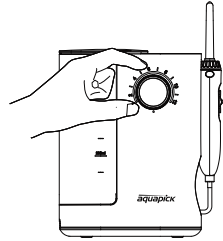
#### 5. Using Oral Irrigator

The product automatically starts working when you lift up the handle. For better experience, when the handle is picked up, the motor will start after 1 second to give you time to put the jet tip into your mouth. Use the handle switch up and down to open and close the water stream.

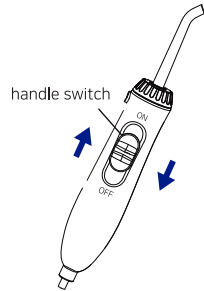
#### 6. Irrigation Guide

- ① Start oral irrigation by pushing the handle switch upwards.
- ② Direct the water stream gently along the teeth.
- ③ Position the water stream at 90-degree angle to the teeth and gums and spray the water stream along the gum line.
- ④ When cleaning the crown, bridge, or between the teeth and the orthodontic bracket, it is more effective to intensely spray at the specific area.

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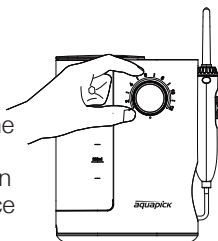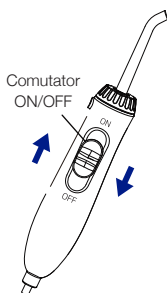
## 7. Pause Control

You can push the handle switch downwards to pause the spraying of water

## 8. When Finished

Attach the handle back to the unit to turn it off.  
Empty any liquid left in the water tank and unplug the unit.

## 10. 3 Minute Timer

The device will automatically switch off after running for 3 minutes. To start again, please put back the handle to the main unit, then pick the handle up again to turn the power back on.

## MAINTENANCE

1. After using the product, invert the water tank until it's dry and store it in a cool and dry place.
2. To wipe clean the product, use soft cloth and neutral detergent.
3. When there is something wrong with the product, please consult the seller.
4. Do not disassemble the product by yourself, or the failure caused by this action will not be covered by the warranty.

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## TROUBLE SHOOTING

Problem	Possible Cause	Solution
Water leakage between the jet tip and the handle	Jet tip hasn't been installed properly	Insert the jet tip into the handle until you hear a clicking sound
Device does not turn on	The engine doesn't work	Ensure that the power source is connected, or contact the retailer for repairs
Can't find the power button	It has no power button	Pick the handle up to turn on the unit and put the handle back on the body to turn the unit off
Device suddenly turns off	It may be turned off because of the 3 minute timer function	Place the handle back and pick it up again
Water leakage from bottom	Some water has been in jet tip holder	Try to prevent water from going into the jet tip holder
	Water tank is not installed properly	Re-install the water tank properly and firmly press it downwards

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Problem	Possible Cause	Solution
Water pressure is weak	Water tank isn't installed properly	Re-install the water tank, and firmly press water tank downwards
The motor works but water doesn't come out from the jet tip	It takes time for water to move from the water tank to the jet tip	Please wait for 1~2 minutes. If water still does not come out, please place the handle back on the device and pick it up again.
	The water pressure is set too low	Please adjust the water pressure to a higher level.
It doesn't work even if I lift the handle.	It could be a one-second delay.	Please wait for 1 second. A second is the time it takes to take the handle to the mouth.

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## IMPORTANT SAFEGUARDS

### ※ Read these instructions before using.

- ① If you are receiving oral treatment, ask your doctor about the use of this product.  
(In the case of side effects during use, stop at once and seek medical advise).
- ② Use product only for its intended use as described in this manual.  
(Do not use for the washing/cleaning of eyes, nose or ears).
- ③ Both tap water and ionic water can be used for the product.  
**※ Do not use water over 40 degrees celcius and be careful not to overfill.**
- ④ Do not use the product for over 7 minutes straight.  
After operating over 7 minutes straight, stop the product for 10 minutes for the motor to cool down before using it again.  
(The motor may overheat and shorten the life of the product.)
- ⑤ Do not place or store the product where temperature is high and where there may be heat.
- ⑥ Do not drop or give external shock to the product.
- ⑦ There may be noise and shaking of the product during use and they come from the motor that pulsates water stream strongly at 1,650 pulses per minute. This is totally normal.

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It is a free voltage product.  
(100~230V, 50/60Hz)

Do not use salt or gargle  
other than water.  
(However, it is recommended  
to use Aquapick Toothfoam)

Do not disassemble,  
reassemble or repair  
the product.


Ensure that water  
does not enter the  
inner part of the body.

Do not use while  
taking a shower.

Close supervision is  
necessary when this  
product is used by  
children.

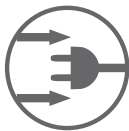Always unplug this  
product immediately  
after using.

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