aquapick[®]

SONIC ELECTRIC TOOTHBRUSH

User Guide

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AQ-101

Thank you very much for using the Aquapick Sonic Electric Toothbrush.

Important safety notes

- When using this product for the first time, due to the high-speed swing of the bristles, if there is a feeling of irritation or tingling, please don't worry. This is a normal phenomenon, and this feeling will gradually decay and slowly disappear.
- If it has not been used for a long time (the battery is out of power due to standby power consumption), the charging indicator is not displayed when charging the toothbrush. Please wait 3-5 minutes and it will resume normal charging.
- If the body of the toothbrush feels slightly hot when charging, this is normal. If it feels very hot, please turn off the power immediately. Then, contact Aquapick customer service.
- This product can only use the special brush head of Aquapick sonic toothbrush .
- This product is only used to clean teeth, gums and tongue. Do not use it for other purposes.
- This product is for personal use only and should not be used in a dental clinic or institution, or for multiple patients at the same time.
- This product is for indoor use only. Do not place the toothbrush in direct sunlight. It should be placed in a dry or ventilated place at normal temperature.
- This product is not intended for the disabled, mentally ill or inexperienced people (including children (1-12 years old)) unless they are used under the supervision of their guardian.
- Please turn off the power when the charging stand is not in use or before cleaning.
- Regularly clean the power plug to prevent the accumulation of dust and dirt, which will reduce the charging effect.

Orally sensitive person

- ❶ People who have just undergone oral periodontal surgery or are being treated.
- ❷ People whom with other symptoms in the mouth, including those with weak dentures or persistent bleeding in the gums.
- ❸ Other people with sensitive oral vulnerability and so on.

Please consult a professional dentist before using this product.

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

Product introduction

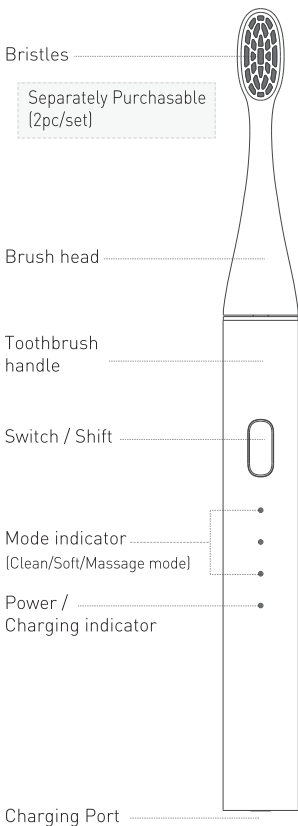
- Start-up: Press the switch button only once and the toothbrush works.
- Switching mode: Within 5 seconds after boot-up, switch mode can be switched by pressing the switch button.

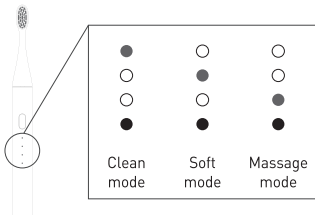
- Mode Automatic Memory: The system records any mode for more than 5 seconds. The memory mode is activated by default on the next boot.
- Shutdown: After 5 seconds, press the switch button; or after 2 minutes, the device shuts down automatically.

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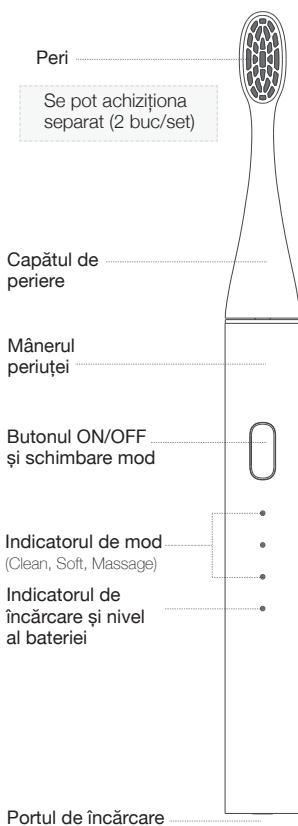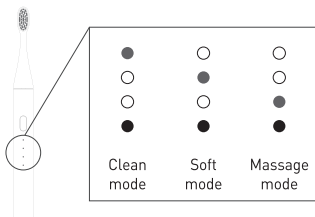
- Charging:

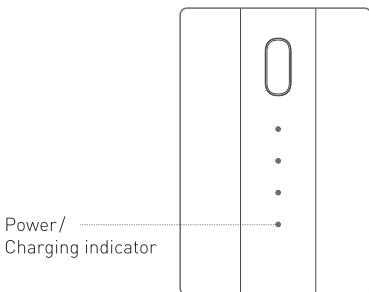When the charging port is in contact with the dedicated charging base, the handle charging indicator lights up, indicating that charging is normal.

Charging Process

- Mode:
 - 0%~10% - indicator is steady red
 - 10%~100% - indicator is steady yellow

The charging is completed and the indicator is steady green.

First Time Using

Please turn on the toothbrush. If there is any power, please use it and recharge when it runs out.


Use Process

If the battery level is lower than 10%, the indicator will light up in red to remind you to charge.

Precautions

- ❶ Before charging, please fully remove the water from the charging port and charger. It could cause corrosion of the product.
- ❷ The battery can not be replaced by the user. It can only be replaced by the Aqapick authorized customer center. Non-professionals are not allowed to disassemble the product by themselves. Otherwise, the user will be deemed to be automatically de-protected.

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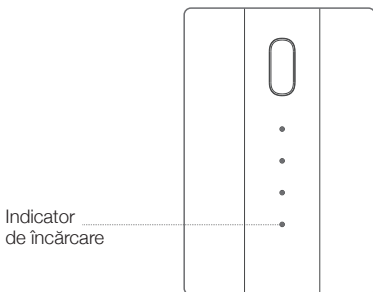
Dental recommended brushing method

Brushing Time


- 1 Dental specialist (dentist) recommendation: brush your teeth with an electric toothbrush, which at least takes 2 minutes to brush each time.
- 2 When you use it, you will feel a short pause every 30 seconds, which reminding you that you should clean the four areas of the mouth to achieve a scientific cleaning effect.

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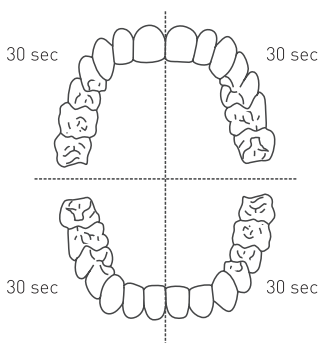
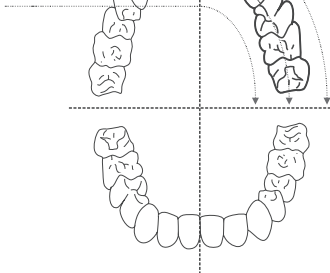
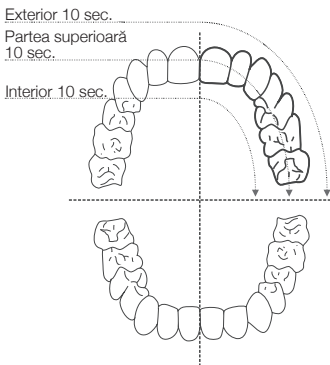
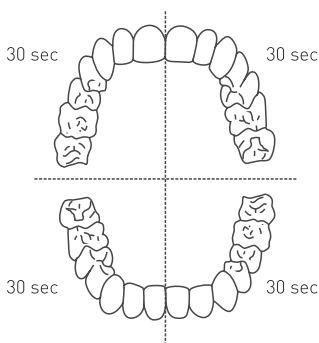
Use Steps

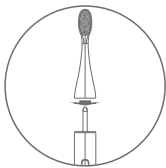Carefully insert the brush head into the swinging shaft, and leave a certain gap between the brush head and the main body to reserve the working space for the brush head to work normally.

To prevent the toothpaste from splashing, take a small amount of toothpaste on the brush before using (Aquapick sonic electric toothbrush saves toothpaste more than manual toothbrush).

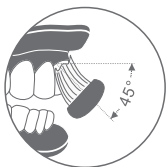Move the bristles to the gum line at 45° and slowly move the toothbrush.


Hold the handle of the toothbrush gently and move it slowly in a circular motion. Do not press it too hard.

Please rinse your mouth with water after brushing.

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
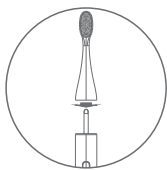
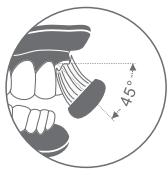
Cleaning and Maintenance

Cleaning brush head

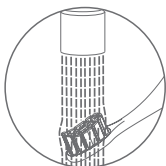
Wash the brush and bristles after each use and allow to air dry. Do not pull the bristles, which will deform them and cause them to fall off.

Cleaning the fuselage

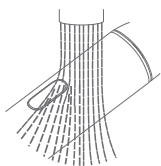
Rinse the handle surface with water and scrub with a soft cloth and air dry. Do not press the rubber seal on the metal shaft with sharp objects as it may cause damage.

Cleaning charger

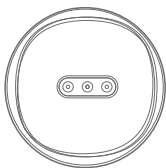
Wipe the dirty cloth on the charger with a soft cloth.
* Always unplug the power cord from the outlet before cleaning the charger. Wipe the surface of the charger with a damp cloth and allow it to dry thoroughly before putting it into service.

Replace the brush head

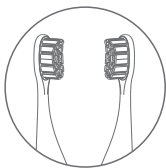
The dentist recommends replacing the brush head every 1-2 months for best results.

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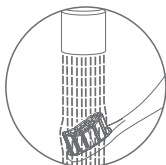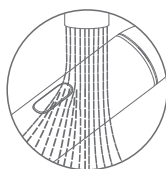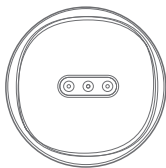
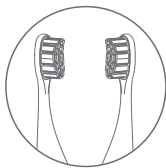
Product Specifications

- Product: Aquapick Sonic Electric Toothbrush
- Model: AQ-101
- Rated voltage/power: 3.7V / 1.85W
- Charging: Aquapick dedicated USB charger
- Input: 5V / 300mAh
- Working hours: 30 Days (*)
- Battery capacity: 800 mAh (**)
- Charging time: 3 Hours (***)

* Calculated according to 2 minutes/time, 2 times/day working condition.

** The capacity of the battery has a tolerance range of ± 50 mAh.

*** Recommended to use the mobile phone industry 5V/1A (or 5V/2A) standard USB charger.

FAQ

The following FAQs will help you use the warranty service effectively. In order to save your time and unnecessary expenses, please do the following simple check before contacting professional after-sales service.

❶ The swing strength is weakened.

Check point:

- Check if the toothbrush is strong on the teeth.
- Check if the gap between the brush head and the body is too large or too tight.

Solution:

- Gently place the brush head on your teeth.
- Remove the brush head from the arm and re-seat it on the output shaft. The tooth head should be kept at a certain distance from the handle.

❷ The toothbrush can only work for a short time.

Check point:

- Check if the charging base is in good contact with the charging interface at the bottom of the main body, resulting in abnormal charging and whether the handle charging indicator lights up red.
- Check for insufficient power.

Solution:

- The charging time of the cleaning base or the charging interface at the bottom of the main body is not enough.
- The battery has reached the end of its life.

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