# **aquapick**<sup>®</sup>

CORDLESS ORAL IRRIGATOR

## **User Guide**

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| **AQ-207**

AQ-207 is a rechargeable Oral Irrigator.

Please notice the following issues before use.

## DANGERS

- 1 If you do not understand the operating procedures and safety instructions, do not use this product to avoid injury caused by improper operation.
- 2 Please do not operate the product in any way other than this manual.
- 3 Do not disassemble the product or change the parts of the product to avoid damage to the product or high voltage electric shock. These acts will make the warranty invalid.

## WARNINGS

- 1 This product is only used for teeth cleaning and gum massage. Please do not use it for other purposes.
- 2 Do not intrude the product into the water to avoid failure, fire, electric shock, explosion and other hazards.
- 3 If the USB cable falls into water, dry it before use.
- 4 Do not plug in or pull out the power plug with wet hands.
- 5 Do not damage or modify the power cord. Do not stretch, bend or twist the power cord. Do not put weight on the power cord.
- 6 Please unplug the USB cable before maintenance.
- 7 Do not place the product in a hot or humid place.
- 8 Please consult the dentist for further use if there is a massive bleeding in the gums after use or if there is still bleeding or other discomfort after one week of use.

## ATTENTIONS

- 1 This product is not suitable for children, infants and people who do not know how to operate the product unless they are supervised or instructed by their guardian. Please keep the product out the reach of children and infants. It cannot be used as a toy.
- 2 Please contact the manufacturer or the dealers if you need a repair service.
- 3 If the toothpaste contains peroxide or bicarbonate (a common ingredient in whitening toothpaste), please clean thoroughly after each use to prevent plastic cracking.
- 4 If you have had oral surgery in the last two months, consult your doctor before using the product.
- 5 This product meets the safety standard of battery devices. If you are using a pacemaker or other implant, consult your doctor before using this product.

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
## CONFIGURATION

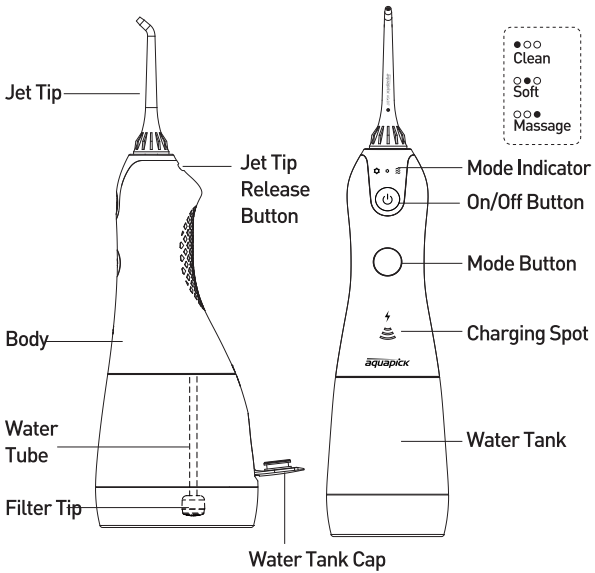
### Note:


There might be some liquid remaining inside the main appliance before use (tank, water suction tube, etc). This is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

## HOW TO CHARGE

Below devices can be used to charge this product. (Fig. 1)

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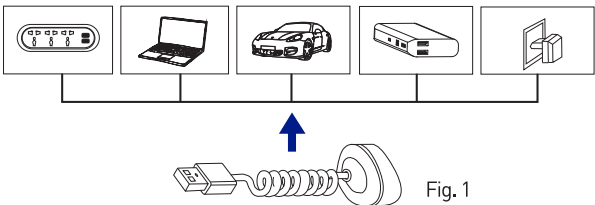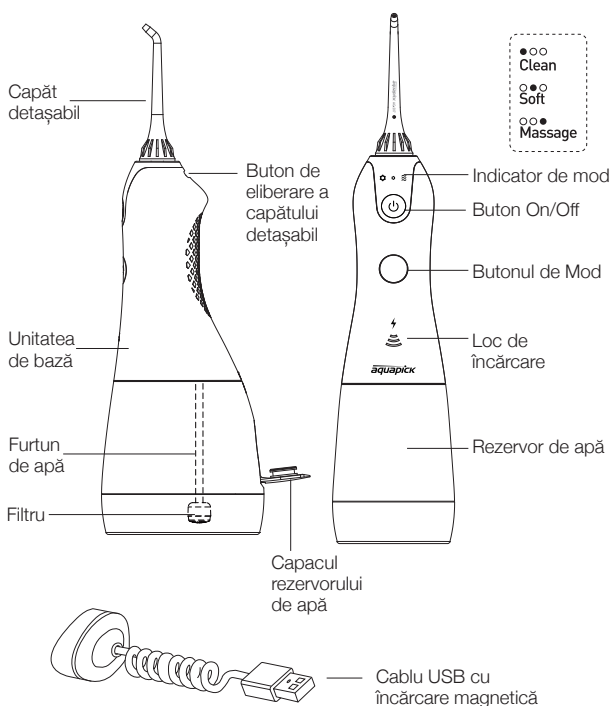
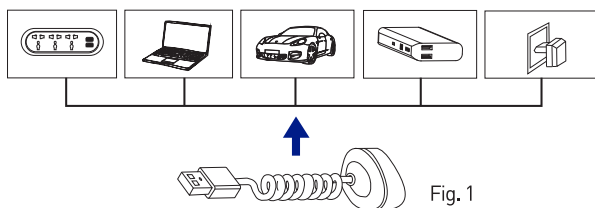
**Note:**

If you place the magnetic part of the USB cable on the charging spot, it will stick to the main body. (Fig. 2)

- When the battery is fully charged, the mode indicator stops blinking and remains white to indicate that it is fully charged.

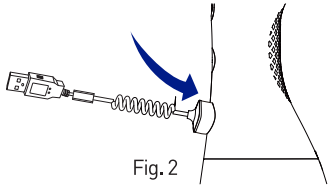
- It takes 12 hours to fully charge the battery.



## USE THE ORAL IRRIGATOR

- 1 Connect the jet tip onto the body. When replacing the jet tip, turn the product off and press jet tip release button to easily remove the current jet tip. (Fig. 1)


- 2 How to put on the water tank (Fig. 2)  
Insert and match the body to the water tank; rotate the water tank by 15° and it will lock onto the body.

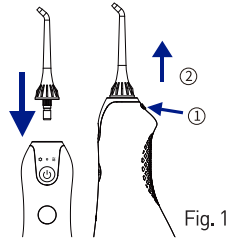- 3 How to remove the water tank (Fig. 3)  
Rotate the water tank by 15° and dislocate it from the body; hold the water tank still and pull up the body.

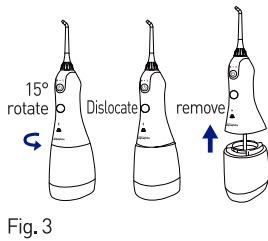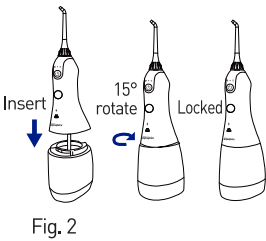- 4 Separate the water tank from the body, fill the tank with water that does not exceed 40°C, and connect it back onto the body. (Fig. 4) Or open the water tank cap to fill the tank with water that does not exceed 40°C, and then close the cap.(Fig. 5)

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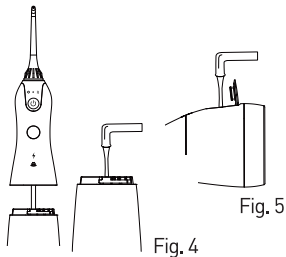Fig. 4

## 5 How to operate

Press the mode button to switch the mode in order of Clean → Soft → Massage. After selecting the mode you like, press the On/Off button to start working. Press On/Off button again to stop working. You can change the operating mode while oral irrigating.

Press the mode button to choose the mode you like.

- Clean Mode: Normal cleaning
- Soft Mode: Gentle cleaning for sensitive gums
- Massage Mode: Cleaning and massaging the gums with pulsating water flow

6 Keep the tip of the jet tip straight toward your teeth.

7 Slightly keep your mouth open so that the water can flow out.

8 Press On/Off button to start irrigating.

1) Make sure the jet tip forms a 90° angle with your teeth and gums.

(Fig. 5)

2) Spray water gently along the teeth.

3) Wash your teeth along the gum line so that the jet tip and gums are at a 90 degree angle.(Fig. 6)

4) When cleaning crowns, bridges, or between teeth and orthodontic brackets, it is better to spray intensively on the affected areas.

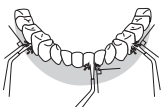
Fig. 5

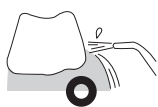
Fig. 6

9 This product has a 2-minute timer function. After 2 minutes of use, the product turns itself off to prevent the motor from overheating.

## MAINTENANCE

1. Turn the product off by pressing On/Off button.

2. Open the water tank cap and dispose of the water in the water tank. (Fig. 1)

3. Press On/Off button to turn on the product and drain the water off completely. (Fig. 2)

4. Press the On/Off button again to turn the product off and close the water tank cap.

5. Wipe the product with dry cloth. (Fig. 3)

6. Press the jet tip release button to remove the jet tip.

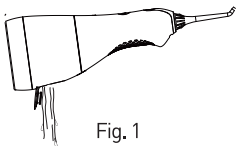
Fig. 1

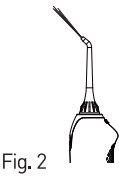
Fig. 2

Fig. 3

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Fig. 3

## NOTE

- To prevent bacteria from spreading, please check whether there is any water remaining inside the product.
- If you are not going to use it for a long time, be sure to wipe and dry the product before storing it.
- Use only water or neutral cleaning agent to clean the product.
- Never use caustic or abrasive agents (e.g: vinegar or lime scale removers), since they may damage the device.
- Do not wash the product with hot water over 40°C.
- Keep the product away from high temperature area and avoid direct sunlight.
- For hygienic use, it is recommended to replace the jet tip every 6 months. (Purchasable separately)
- Do not bend, pull, or twist the jet tip.

## TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
It does not work	You've just purchased the product or it has been left unused for more than 3 months	Charge the unit for about 4 hours
It only works for a few minutes even after charging	The battery has reached the end of its life cycle	Contact repair service center
	Charging time was insufficient	Charge the unit for about 4 hours
Water pressure is very low	There is something stuck in the jet tip	Replace the jet tip
	The jet tip is deformed	
	There is something stuck in the filter	Clean the filter or purchase the filter separately
Water does not come out	The water tank is empty	Fill up the water tank
	The product is tilted too much while in use	Use the product in a vertical state
The mode indicator does not light up when charging	Power is not supplied	Connect the USB cable properly
	The USB cable is defective	Contact repair service center

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