aquapick[®]

PORTABLE ORAL IRRIGATOR

User Guide

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IMPORTANT SAFE GUARDS

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

READ ALL INSTRUCTIONS BEFORE USING.

1. DANGERS

- 1 Toddlers, infants, and people who are unable to operate should not be allowed to use this appliance.
- 2 Never alter the appliance. Also, do not try to disassemble or repair it.
- 3 Do not place or store the product where it can fall or be pulled into tub or sink.
- 4 Do not use the product while bathing.

2. WARNINGS

- 1 Never use the adaptor if it has a damaged cord or plug.
- 2 Keep the cord away from heated surfaces.
- 3 Please do not use it if you have lost the jet tip.
- 4 Do not fill the tank with water warmer than 40°C.
- 5 Do not submerge the device in water.
- 6 Do not use for any purpose other than oral cleaning.
- 7 Keep the jet tip away from small children. The jet tip is small enough to be swallowed and can be lodged in the throat.
- 8 People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- 9 Do not press the jet tip too firmly against your teeth or gums. Doing so may injure your gums.
- 10 This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities or lack of experience and knowledge unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety.
- 11 The appliance is only to be used with the power supply unit provided with the appliance.
- 12 Do not attempt to remove the battery.
- 13 The battery must be removed from the appliance before it is scrapped.

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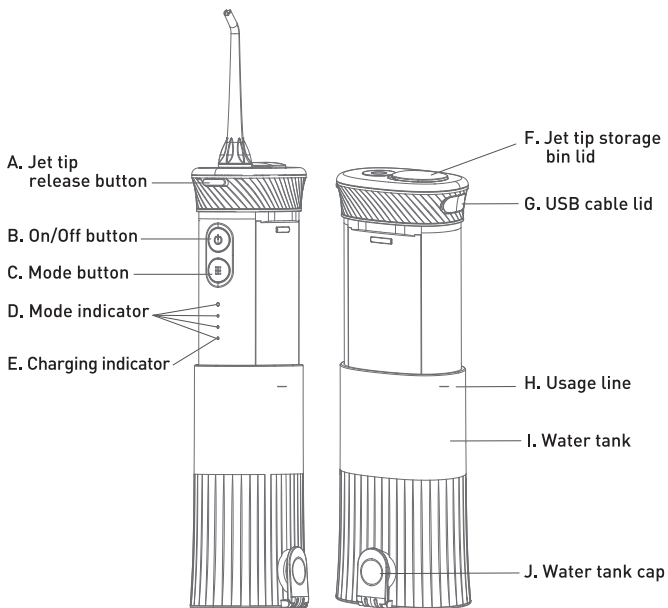
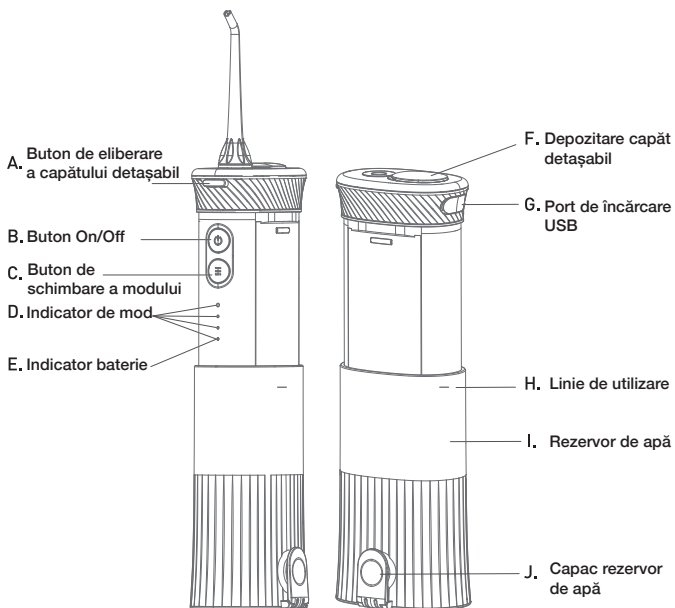
HOW TO USE

Charging the device (Fig. 1/2)

1. Take the charge socket cap out and insert the charging cable securely into the charge socket (make sure no liquid remains inside the socket before inserting the charge socket).
2. Insert the charging adaptor securely into the power outlet.
3. The charging indicator blinks when the device is being charged. It will stop blinking and the light will stay on once the device is fully charged.
4. Put the charge socket cap back into the charge socket after charging is finished.

Note :

If the charging indicator keeps blinking during or after irrigating, it means that the battery needs to be charged. It takes about 4 hours to fully charge the device.

Do not turn on the device before placing the jet tip in your mouth.

1. Open the water tank cap. (Fig.3)
2. Pull out the main unit up to the usage line. (Fig.4)

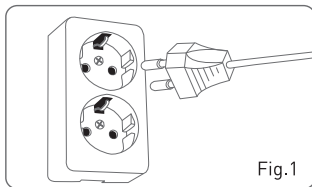
Fig.1

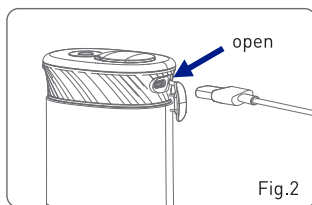
Fig.2

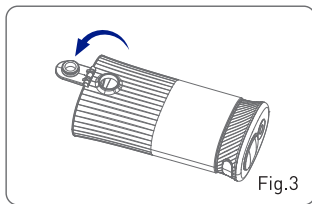
Fig.3

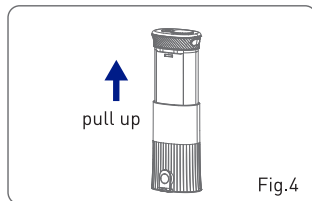
Fig.4

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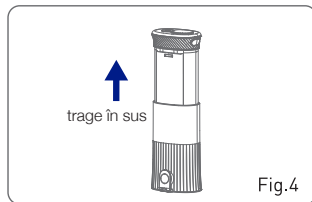
Fig.4

3. Attach the jet tip onto the unit. (Fig.5)

4. Fill the water tank with water which temperature not exceed 40°C, then close the water tank cap. (Fig.6)

5. Press "Mode" button to select the mode you desire. During irrigation, you can change the operating mode.

- Clean: Normal cleaning
- Soft: Gentle cleaning for sensitive gums
- Massage: Cleaning and massaging the gums
- Custom: Freely select desired water pressure to clean and massage the gums

* How to change into Custom mode

After selecting Custom mode by pressing Mode button, press Power button. After that, press and hold Mode button for 3 seconds until the mode indicator blinks. Each press increases the water pressure from 1 to 7 levels, giving you the freedom to choose the water pressure you want.

6. Hold the unit in an upright position with nozzle tip point toward your teeth.

7. Keep your mouth opened slightly so that the water can flow out. (Fig.7)

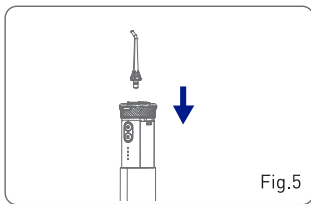
Fig.5

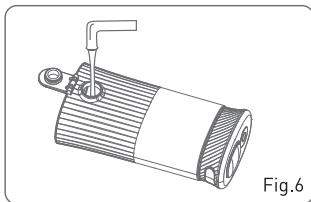
Fig.6

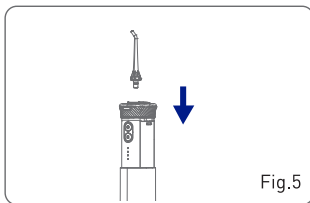
Fig.7

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
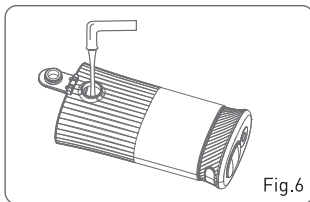
8. Press "On/Off" key to start irrigating. (Fig. 7)


- 1) Directing the water jet flow at a 90-degree angle to the teeth and gums.
- 2) Directing the water jet slowly along the teeth.
- 3) Direct the water stream along the gum line with the jet tip at a 90-degree angle to the gum tissue.
- 4) Clean between the teeth and around any braces, crowns, or bridge work thoroughly
- 5) Press "On/Off" button and then press the jet tip release button to replace other jet tip when necessary.

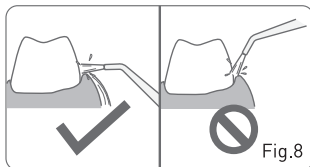- 6) Once irrigating is finished, press "On/Off" button again to turn off the device.

After using the appliance (Fig.9)


1. Open the water tank cap and dispose of the water in the tank.
2. Press "On/Off" button to turn on the unit, and drain the water inside the appliance.
3. Press "On/Off" button again to turn off the unit, and close the water tank cap.
4. Wipe the device with dry cloth.

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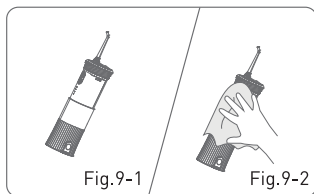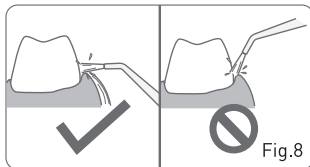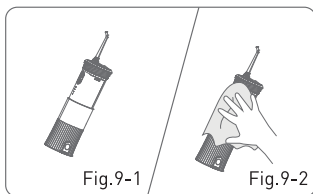
MAINTENANCE

Note:

To prevent bacteria propagating, make sure there is no remaining water inside the unit. If you're not going to use this appliance for long time, be sure wipe and dry the appliance ASAP before store it.

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 50°C. Keep the device away from high temperature area and avoid direct sunlight.

Main unit

1. Remove the dirt by wiping using a cloth with liquid soap and wipe off using a dry cloth.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.

Water tank cap

Wipe off any grime and dirt with a cloth. Water suction tube

Jet tip

1. Wash with water and wipe with soft cloth.
2. Replace the jet tip every 3 months.

Water suction tube

1. Wash with water and wipe with soft cloth.
2. Do not bend, pull or twist the water suction tube.

Water tank

1. Wash with water and wipe with soft cloth.
2. Dry the tank when not using for a long period.

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TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
It does not work	You've just purchased the product or it has been left unused for more than 3 months.	Charge the unit for at least 4 hours.
It only works for a few minutes even after charging	The battery has reached the end of its lift cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for at least 4 hours.
Water pressure is very low	There is dirt stuck in the jet tip	Replace the jet tip
	The jet tip is deformed	
	There is dirt stuck in the filter	Clean the filter
Water does not come out	The tank is empty	Fill up tank
	The appliance is tilted too much while in use	Use with the appliance in a vertical state
	You've just purchased the product or it has been left unused for a period.	Pull out the jet tip, drop a few drops of water in the water hole and operate the product.

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



DISPOSAL

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below. During below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

1. Remove screws lid and screws on the bottom the handle (Fig. 10)
2. Use a sharp stick to pry the top cover open (Fig. 11).
3. Open the shell and take out the core components (Fig. 12).
4. Cut off the battery line and pull out the battery(Fig. 13).

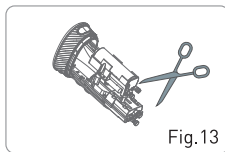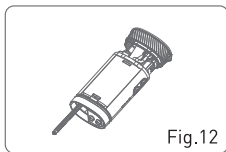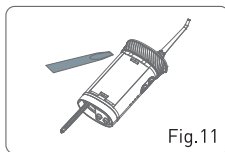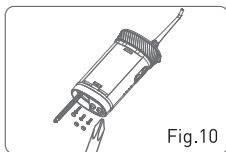If you have recycling questions, please contact your local waste management office.

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